


medilogAR


Efficient and reliable state-of-the-art Holter

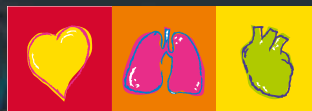

SCHILLER
The Art of Diagnostics

medilogAR

- Robust, shock- and splash-proof as well as easy to clean
- Flexible dual-battery concept for more than 14 days recording duration
- Instantaneous atrial fibrillation/atrial flutter detection based on true P-wave analysis^{1,2}

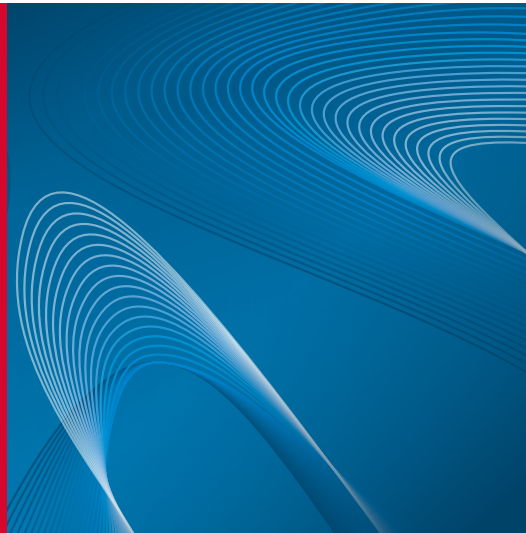

IDEALLY SUITED FOR CLINICAL ROUTINE

The shock- and splash-proof design makes the medilogAR Holter highly suitable for the daily clinical routine. The Holter and its carrying case are robust and easy to clean.

TWO WEEKS OF POWER

Thanks to the dual-battery concept, the patients can be screened for more than 14 days without having to see the doctor to change the batteries.



Eco-friendly dual-battery concept (illustration)

ECO-FRIENDLY CHOICES:

The high-quality internal battery is recharged via USB. By adjusting the signal resolution, the internal battery can run for 24/72 hours or more than 7 days.

AAA SINGLE USE BATTERY:

For fastest turnaround times in high throughput clinics, use easily available standard AAA batteries.

MORE THAN 14 DAYS:

For maximum recording duration, combine a fully charged internal battery with an AAA battery to get more than 14 days with automated AF detection^{1,2}.

VARIOUS CARRYING OPTIONS

medilogAR is so flat that it can easily be carried in a trouser pocket. With the carrying case, it can be worn on a belt or with a neck belt directly on the body.

MADE TO SUPPRESS ARTIFACTS

The medilogAR Holter recorder features a superior resting ECG grade 32,000Hz sampling rate on 3 channels. This is used for true P-wave analysis^{1,2}, sophisticated artifact suppression and motion detection, enabling a speedy evaluation.

INSTANTANEOUS ATRIAL FIBRILLATION / AFL DETECTION:

The medilogAR scans for atrial fibrillation and atrial flutter by directly looking for P-waves^{1,2}, just as a cardiologist would.

¹ Available with medilogAR Professional

² Available with medilogAR Enterprise



RESPIRATION EPISODES:

Thanks to ECG-derived respiration recording², the medilogAR is able to screen for potential respiratory episodes during sleep. The optional SpO₂ sensor², connected via Bluetooth, allows for additional respiratory information.



SpO₂ sensor, attached for the night

MUCH MORE THAN ARRHYTHMIA SCREENING

Instantaneously locating Atrial Fibrillation^{1,2}, analysis of respiratory episodes², and high-resolution HRV are just some of the exceptional features of the medilogAR. The analysis is carried out using the renowned medilog DARWIN2 software.

1 Available with medilogAR Professional

2 Available with medilogAR Enterprise



Americas
SCHILLER Americas Inc.
 Doral, Florida 33172

North America:
 Phone +1 786 845 06 20
 Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com

Latin America & Caribbean
 Phone +1 305 591 11 21
 Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com



Asia
SCHILLER Asia-Pacific / Malaysia
 52200 Kuala Lumpur, Malaysia
 Phone +603 6272 3033
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
 A-4040 Linz
 Phone +43 732 70 99 0
sales@schiller.at
www.schiller.at



China
Alfred Schiller (Beijing) Medical Technology Co., Ltd.
 100102 Beijing, China
 Phone +86 10 84565453
info@schillerchina.com
www.schillermedical.cn



Croatia
SCHILLER d.o.o.
 10000 Zagreb
 Phone +385 1 309 66 59
info@schillerzg.hr
www.schiller.hr



France
SCHILLER Médical
 F-67160 Wissembourg
 Phone +33 3 88 63 36 00
info@schiller.fr
www.schiller.fr



France (distribution France)
SCHILLER France S.A.S.
 F-77608 Bussy St Georges
 Phone +33 1 64 66 50 00
contact@schillerfrance.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
 D-85622 Feldkirchen b. München
 Phone +49 89 62 99 81 0
info@schillermed.de
www.schillermed.de



India
SCHILLER Healthcare India Pvt. Ltd.
 Mumbai – 400 059, India
 Phone +91 22 6152 3333 / 2920 9141
sales@schillerindia.com
www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
 PL-02-729 Warszawa
 Phone +48 22 647 35 90
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
 119049 Moscow, Russia
 Phone +7 (495) 970 11 33
mail@schiller.ru
www.schiller.ru



Serbia
SCHILLER d.o.o.
 11010 Beograd
 Phone +381 11 39 79 508
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
 2310 Slovenska Bistrica
 Phone +386 2 843 00 56
info@schiller.si
www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
 ES-28232 Las Rozas/Madrid
 Phone +34 91 713 01 76
schiller@schiller.es
www.schiller.es



Switzerland
SCHILLER-Reomed AG
 CH-8912 Obfelden
 Phone +41 44 744 30 00
sales@schiller-reomed.ch
www.schiller-reomed.ch



Türkiye
SCHILLER TÜRKIYE
 Okmeydani-Sisli – Istanbul
 Phone +90 212 210 8681 (pbx)
 Fax +90 212 210 8684
info@schiller.com.tr
www.schiller-turkiye.com



United Kingdom
SCHILLER UK Ltd.
 Anstruther, Fife KY10 3H
 Phone +44 133 331 2150
sales@schilleruk.com
www.schilleruk.com

SCHILLER reserves the right to make changes in design, specifications and features shown, or discontinue the product described at any time without notice or obligation. Images may be representative. Some products may not be regulatory cleared or released for sale in all markets. Contact your SCHILLER representative for the most current information and availability.



SCHILLER
 The Art of Diagnostics